

Commander Echo Inclinometry for Spine and Extremity Range of Motion

Wireless Commander Echo inclinometers provide more features and functionality than any non-computerized inclinometer system. It makes spine and extremity range of motion measurements faster and easier.

Advanced Features

- Measure true range of motion plus movement at the upper and lower extremes of the spine segment without resetting the inclinometer placement
- Automatically collects range of motion end points without pressing a button in dynamic testing mode
- Perform a complete spine Range of Motion (ROM) exam without pausing to write down results
- Quantifies total range of motion and motion of each inclinometer
- Automatically subtracts extraneous movement for true spine ROM measurements
- Static test mode easily simulates most goniometer ROM techniques
- Wireless console is independent of the testing device allowing the tester to review data during active testing without influencing the patient's level of effort



*Console required to operate device. Console and device sold separately.

DISTRIBUTED BY



190, Boul. St-Elzéar Ouest
Laval (Québec) H7L 3N3
T 450 669.8985 | 1 800 361.0877
F 450 669.9532 | 1 888 935.7001
www.medicaltronik.ca



JTECH MEDICAL
MEASURING HUMAN PERFORMANCE
© 2012 JTECH MEDICAL MM317_B

Technical Specifications


- Number of tests stored on Console:20
- Number of repetitions per test: 1-8
- Audible cues: yes
- Show independent angle readings yes
- Units of measure degrees
- Allows for specific test re-testing: Last rep only while still testing
- Print options: Yes*

* With purchase of Downloader Echo software

Included Accessories:

Alignment Rails (8MH068) – Magnetically attach to the Inclinerometers for easy visual confirmation of alignment.

VELCRO® Straps (45, 46, 47) – Can be used for dynamic range of motion inclinometry testing. Comes in a set containing Small, Medium, and Large straps.

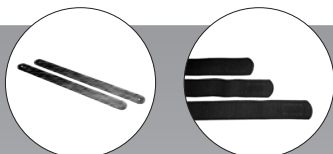

Product Highlights

- Console stores data for 20 tests with up to 8 repetitions each
- Measure lordosis, kyphosis and joint ankylosis
- Measure agonist/antagonist movements, like flexion and extension, without stopping or repositioning
- Achieves accuracy and repeatability of +/- 2 degrees
- Measure ROM using dynamic, auto-rep or static testing

COMMANDER ECHO Downloader™

Optional Commander Echo Downloader (CM311) software enables you to create quick but meaningful reports from the

data you collected on your Commander Echo Console with the click of a button. The easy-to-use Commander Echo Downloader software interfaces directly with your Commander Echo Console eliminating dictation and transcription errors.



Other Commander Echo Devices*



*Some devices shown with accessories attached.

DISTRIBUTED BY



190, Boul. St-Elzéar Ouest
Laval (Québec) H7L 3N3
T 450 669.8985 | 1 800 361.0877
F 450 669.9532 | 1 888 935.7001
www.medicaltronik.ca



JTECH MEDICAL
MEASURING HUMAN PERFORMANCE

© 2012 JTECH MEDICAL MM317_B