



A Safer Way for Both Caregiver and Client

The Hera highlights the best in European design and attention to detail. Hera is a sophisticated yet easy to use bathing system that promotes bathing and wellness. Hera, known by its compactness, is freestanding and easy to install. Hera is constructed utilizing the utmost in European fiberglass craftsmanship.

No compromise on infection control standards, the dual hydromassage jets utilize independent, self-contained, whirlpool jets. There is no plumbing or tubing running between jets and no central whirlpool pump.

Contact Waverley Glen the specialists in health care facility design, and moving and handling to learn more about our solutions.

Product Highlights

- "Easy Entry" electronically locking and opening side door, opens a full 180°.
- Accessible by total lifts, sit/stand lifts, ceiling lifts, chair lifts or ambulatory clients.



- Fully automated panel, one touch disinfection system can be added to prevent cross infection hazard.
- Optional electronic dispensing of bath additive and shampoo, which is located in a lockable storage compartment.




- One touch activation of pre-selected water fill and shut off.



Lift • Position • Transfer • Solutions

A Prism Medical Company

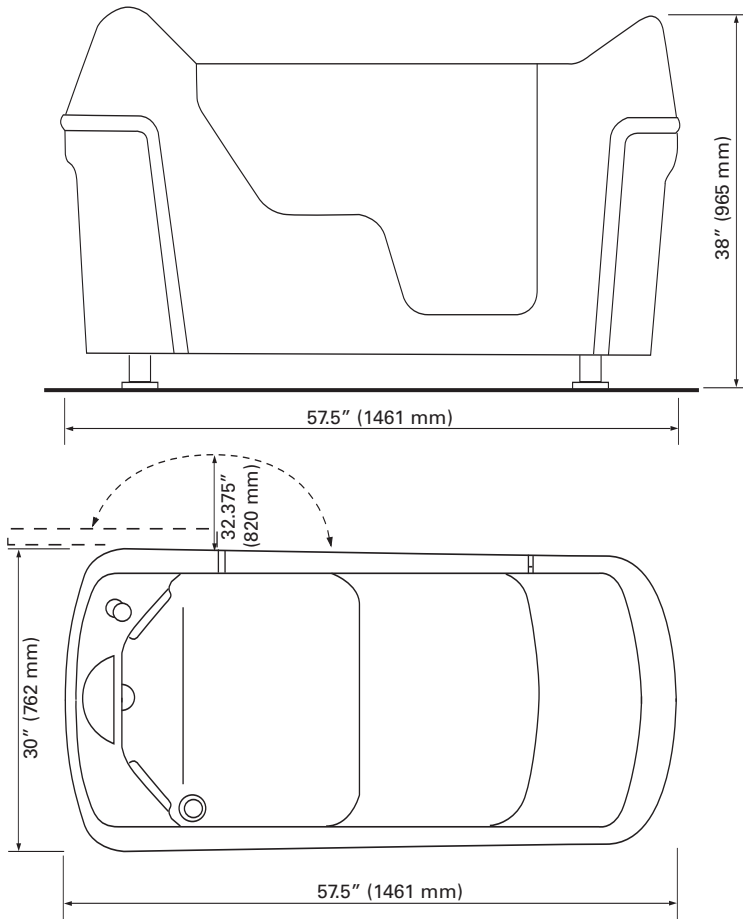

Technical Specifications

- Maximum Width:** 30" (760 mm)
- Maximum Length:** 57.5" (1461 mm)
- Maximum Height:** 38" (965 mm)
- Water Connections:** 3/4" drain Ø 50mm
- Drain Connections:** Flexible Hose
- Capacity of Tub:** Approximately 47.5G (180L)
- Electrical Requirements:** 115VAC, 15 amp GFI
- Safety Class:** I IP 24

DISTRIBUTED BY



MTM 190, boul. St-Elzéar O., Laval (Qc) H7L 3N3
 Tel: 450 669-8985 □ 1 800 361-0877
 MÉDICAL TRONIK Fax: 450 669-9532 □ 1 888 935-7001
www.medicaltronik.ca

www.waverleyglen.com | 1.800.265.0677

Product Highlights

- Manually adjusted thermostatic mixing valves for filling and showering.
- Digital readout of tub and shower temperature.



- Optional safety grab handles.
- Ergonomically designed comfort seat and sloped foot well for additional leg and foot comfort.



- Maximum strength ultra-durable door components ensure years of trouble free performance.
- Dual backup system for opening the bathtub in absence of power.



Lift • Position • Transfer • Solutions

A Prism Medical Company

116 Rayette Road, Unit 1
 Concord, ON Canada L4K 2G3
T 905.660.0213 F 905.660.5513